


To apply the Grapnel tourniquet the clip is held in one hand with the two knobbed flukes pointing down towards the limb. The stretched rubber tubing is then carried once, twice or thrice round the limb, passed beneath one of the knobbed flukes over the shank and out beneath the other knobbed fluke and let go. The fastening may be rendered still more secure if the rubber, instead of being led out beneath the second fluke, is passed, as in the diagram, further round the shank and brought up between the first fluke and the rubber.

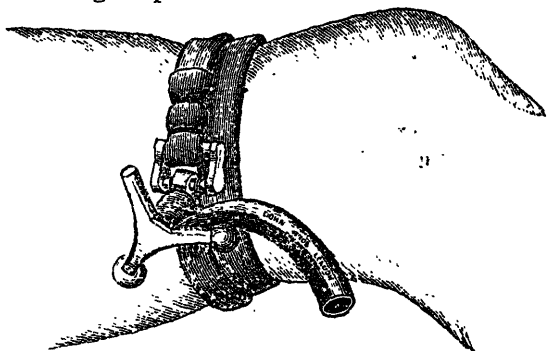

The advantages claimed for these clips are:

1. They are small and light.
2. They are simple and cannot get out of order.
3. They are very quickly applied or removed.
4. The tighter the rubber is applied round the limb the more securely the clip holds it.
5. The clip does not cut the rubber.
6. The clip fastens equally well at any point along the length of the rubber.